



Q-SERIES

**THE
POWER OF
QUICKIE**



Q-Care Range

Advanced engineering. Customised for you.



Q300



Q400



Q500



Q700

ADVANCED ENGINEERING CUSTOMISED FOR YOU



QUICKIE®, SEDEO®, JAY® and WHITMYER® bring you the most comprehensive range of premium, innovative and technically advanced powerchairs, each created individually for you whatever your clinical or lifestyle need.

SUMMARY

THE NEW QUICKIE Q-CARE POWER RANGE: A COMPLETE POWER PRODUCT PORTFOLIO UNDER ONE ROOF

A wide range of bases, drive tractions, seating systems, electronics, controls and accessories, that can be combined to create the perfect powerchair for each individual need.



SERIES

Q700
Q500
Q400
Q300



BASES

Front-Wheel Drive
Mid-Wheel Drive
Rear-Wheel Drive



SEATING SYSTEMS

SEDEO LITE
CAPTAIN SEAT



CONTROLS

VR2 Joystick

SUMMARY

CHOICE WITHOUT COMPROMISE

The new QUICKIE Q-Care power range offers flexibility to create the perfect powerchair for each individual user. Highly flexible configurations, interchangeable parts, quick and easy adjustments.



BEST IN CLASS SEATING

Q-Care offers a range of seating systems to fit clinical needs including Sedeo Lite and Captains Seat. All compatible with the JAY range.

BEST DRIVE BASES CHOICE IN THE MARKET

From cost effective compact bases to high performance outdoor bases, all with the most extensive choice of drive tractions. An impressive range of base levels, and drive tractions that can accommodate all sorts of user preferences, lifestyle needs and pricing/reimbursement requirements.

EASIER TO LEARN

The common design and the modularity create a simplified user experience, that makes it very easy to learn, manage and prescribe.


EASIER TO SERVICE

The modularity and flexibility of the bases, seatings and controls allows an easy adaptation to different users, making it perfect for refurbishments and allowing a great fleet management.

RANGE OVERVIEW

	Q700	Q500	Q400	Q300
Outdoor performance	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Indoor performance	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Suspension	● ● ● ● ● All Wheels	● ● ● ● ● ● Drive wheels, Advanced	● ● ● ● ● ● Drive wheels, Regular	● ● ● ● ● ● Drive wheels, Basic
Motors	4 pole	4 pole	4 pole	2 pole
Speed	10 kph	10 kph	10 kph	10 kph
Batteries	60 Ah / 72 Ah	60 Ah	60 Ah	34 Ah
Max. user weight	136 kg	136 kg	136 kg	136 kg
Drive bases	MWD	MWD / FWD / RWD	MWD / FWD / RWD	MWD / RWD
Seating	CAPTAIN SEAT	SEDEO LITE / CAPTAIN SEAT (MWD only)	SEDEO LITE / CAPTAIN SEAT (MWD Only)	SEDEO LITE / CAPTAIN SEAT (MWD)
Electronics	VR2	VR2	VR2	VR2
Crash tested	Yes	Yes	Yes	Yes

**CREATED INDIVIDUALLY
FOR YOUR CLINICAL AND
LIFESTYLE NEEDS**

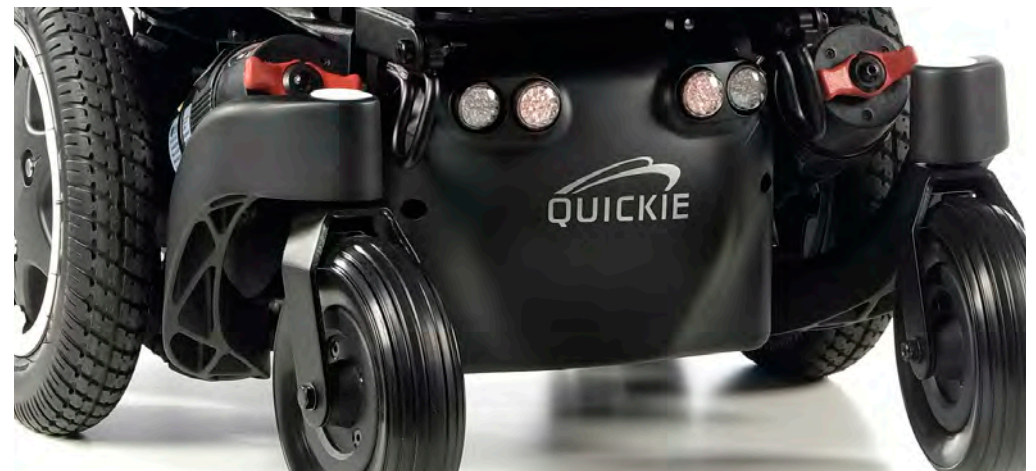

RANGE OVERVIEW

QUICKIE Q-CARE RANGE OVERVIEW

Q700 RANGE: THE ULTIMATE EXPERIENCE

Chairs designed to offer a superb outdoor experience while allowing a great indoor use. Our mid wheel drive base with market leading suspension technologies, powerful motors and long lasting batteries will provide the ultimate outdoor performance.

Q500 RANGE: MAKE YOUR OWN ADVENTURE

Robust outdoor powerchairs with a choice of three drive bases with advanced suspension technologies and powerful motors, for a superior outdoor experience.

Q400 RANGE: ADVANCED OUTDOOR / INDOOR PERFORMANCE

Robust outdoor powerchairs with a choice of three drive bases, with regular motors and suspension for a great outdoor experience with advanced electronics.

Q300 RANGE: ULTRA COMPACT EVERY DAY PERFORMANCE

Entry level chairs with ultra-compact mid or rear wheel drive bases for great indoor agility with the stability and traction for outdoor navigation for community access.





PERFORMANCE BASES

DRIVE BASES OVERVIEW

The new QUICKIE Q-Care range offers the most complete choice of drive bases on the market, available in different series and with many seating systems.

The choice of a drive base affects how the chair performs in different situations and environments. The right choice for each user will depend on personal preferences, experience, clinical, lifestyle and environmental needs.

	MWD	FWD	RWD
Manoeuvrability	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Stability	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Kerb climbing	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Traction	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Drive comfort	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
Ease of driving	● ● ● ● ●	● ● ● ● ● ●	● ● ● ● ● ●
	MWD	FWD	RWD
Q700	✓	—	—
Q500	✓	✓	✓
Q400	✓	✓	✓
Q300	✓	—	✓

PERFORMANCE BASES



MID-WHEEL DRIVE BASES (MWD)

The mid-wheel drive bases have the best indoor performance. Their smallest turning circle provides great manoeuvrability. They are very intuitive and easy to drive, ideal for users with reduced hand control, driving with switches and users with poor spatial awareness.



FRONT- WHEEL DRIVE BASES (FWD)

The front-wheel drive bases have great outdoor performance with good indoor manoeuvrability. Knee angle can be positioned at 90°. Maintains control on steep declines.



REAR-WHEEL DRIVE BASES (RWD)

The rear-wheel drive bases are stable and have good kerb climbing ability. Intuitive drive that provides good control with non-proportional inputs.

SEATING SYSTEMS



SEDEO LITE

Multi-adjustable seat frame with the option of either SEDEO LITE backrest cushions, Jay backrests or sling back upholstery. SEDEO LITE is designed for users with medium to low clinical needs that require adjustable seating system. The backrest, armrests and legrests can be adjusted to fit the user. Includes seat frame with seat rail for additional mounting solutions.



CAPTAIN SEAT

Provides a relaxed sitting position and features a semi-recline, fold-down backrest, adjustable seat size, and convertible seat pan that allows you to use your preferred cushion. Available with power lift or 30°/50° tilt and power centermount legrest.



WHITMYER & JAY SEATING

Compatible with SEDEO LITE seating systems, the JAY family of cushions and backs offers superior clinical seating for pressure management, positioning and comfort. WHITMYER headrests are also compatible with SEDEO LITE seating systems offering the best head support and positioning.



SEATING SYSTEMS

SEATING OVERVIEW

	SEDEO LITE	CAPTAIN SEAT
User clinical need	Medium to Low	Low
Seat width	385 - 510 mm	16 - 20 Inches
Seat depth	410 - 560 mm	16 - 22 Inches
Backrest height	460 - 510 mm	25 inches
Configurability	Medium	Low
Powered recline	85°-120°	—
Manual recline	90°-120°	90°-115°
Tilt angle	30°/50° *	30°/50° *
Lift	300 mm	300 mm
JAY & WHITMYER compatibility	✓	Cushion only
Max. user weight	136 kg	136 kg
Crash tested	Yes, ISO	Yes, ISO

*50° tilt available on Q500 and Q700 only

VR2 JOYSTICK

The QUICKIE Q-Care Series wheelchairs come with the reliable VR2 control. The joystick features easy to use buttons and intuitive mode selection.



VR2 Joystick

ATTENDANT CONTROL

Allows for a carer to assist with driving or power seat function control.



Attendant Control



DESIGN AND APPEARANCE

DISTINGUISHED STYLE

With the Q-Care Range, we've taken style to the next level. With a matte black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and a contemporary design with a variety of colours to choose from.

Fancy a change?

A choice of contemporary colour inlays to perfectly complement its striking matte black bodywork. Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.



*Q700/Q500: 

*Q400: 

*Q300: 



* The number of colours available may differ depending on the seat and base combination.
* Q300 has shroud only, not changeable inlays.



Sunrise Medical Australia Pty Ltd.
6 Healey Circuit
Huntingwood NSW 2148
Phone: (02) 9678 6600
Fax: (02) 9423 6969
Email: enquiries@sunrisemedical.com.au
www.SunriseMedical.com.au

© 2019 Q-Care Range Brochure AU Rev B

For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

Pictures and illustrations may be overseas models and some seating and accessories not available in Australia.

If you are visually impaired, this document can be viewed in PDF format at SunriseMedical.com.au.

